



D-Dimer Fast Test Kit (Immunofluorescence Assay)

IF1006 for Getein1100
IF3006 for Getein1180
IF4006 for Getein1200
IF2006 for Getein1600
IF5006 for Getein1160
IF6006 for Getein208



User Manual

INTENDED USE

D-Dimer Fast Test Kit (Immunofluorescence Assay) is intended for *in vitro* quantitative determination of D-Dimer in human plasma or whole blood samples. The test is used as an aid in the assessment and evaluation of patients suspected of deep-vein thrombosis or pulmonary embolism.

SUMMARY

Deep-vein thrombosis is a common condition, with a lifetime cumulative incidence of 2 to 5 percent. Untreated deep-vein thrombosis can result in pulmonary embolism, a potentially fatal outcome. Anticoagulant therapy reduces both morbidity and mortality from venous thromboembolism, and early diagnosis is therefore important. Accurate diagnosis of deep-vein thrombosis minimizes the risk of thromboembolic complications and averts the exposure of patients without thrombosis to the risks of anticoagulant therapy.

D-Dimer is a marker of endogenous fibrinolysis and should therefore be detectable in patients with deep-vein thrombosis. In recent years, an increasing number of studies have shown the D-Dimer assay has a high negative predictive value and D-Dimer is a sensitive but nonspecific marker of deep-vein thrombosis. Negative D-Dimer can exclude deep-vein thrombosis and pulmonary embolism.

PRINCIPLE

The test uses an anti-human D-Dimer monoclonal antibody conjugated with fluorescence latex and another anti-human D-Dimer monoclonal antibody coated on the test line. After the sample has been applied to the test strip, the fluorescence latex-labelled anti-human D-Dimer monoclonal antibody binds with the D-Dimer in sample and forms a marked antigen-antibody complex. This complex moves to the test card detection zone by capillary action. Then marked antigen-antibody complex is captured on the test line by another anti-human D-Dimer monoclonal antibody. The fluorescence intensity of the test line increases in proportion to the amount of D-Dimer in sample.

Then insert test card into Getein1100/Getein1160/Getein1180 Immunofluorescence Quantitative Analyzer/ Getein208 Hand-

held System (automatically inserted by Getein100/Getein1600 Immunofluorescence Quantitative Analyzer (hereafter referred to as Getein1100, Getein1160, Getein1180, Getein208, Getein1200 and Getein1600)), the concentrations of D-Dimer in sample will be measured and displayed on the screen. The value will be stored in Getein1100/Getein1160/Getein1180/Getein208/Getein1200/Getein1600 and available for downloading. The result can be easily transmitted to the laboratory or hospital information system.

CONTENTS

1. A kit for Getein1100/Getein1160/Getein1180/Getein208 contains:

- Package specifications: 25 tests/box, 10 tests/box
1) D-Dimer test card in a sealed pouch with desiccant
2) Disposable pipet
3) Sample diluent
4) User manual: 1 piece/box
5) SD card: 1 piece/box

2. A kit for Getein1200/Getein1600 contains:

- Package specifications: 2×24 tests/kit, 2×48 tests/kit
1) Sealed cartridge with 24/48 Getein D-Dimer test cards
2) User manual: 1 piece/box
Materials required for Getein1200/Getein1600:

- 1) Sample diluent: 1 bottle/box
2) Box with pipette tips: 96 tips/box
3) Mixing plate: 1 piece/box

3. Sample diluent composition:

Phosphate buffered saline, proteins, detergent, preservative, stabilizer.

4. A test card consists of:

A plastic shell and a reagent strip which is composed of a sample pad, nitrocellulose membrane (one end of the membrane is coated with a fluorescence latex-labelled anti-human D-Dimer monoclonal antibody, the test line is coated with another anti-human D-Dimer monoclonal antibody and the control line is coated with rabbit anti-mouse IgG antibody), absorbent paper and liner.

Note: Do not mix or interchange different batches of kits.

APPLICABLE DEVICE

Getein1100 Immunofluorescence Quantitative Analyzer
Getein1180 Immunofluorescence Quantitative Analyzer
Getein1200 Immunofluorescence Quantitative Analyzer
Getein1600 Immunofluorescence Quantitative Analyzer
Getein1160 Immunofluorescence Quantitative Analyzer
Getein208 Hand-held Integrated System

STORAGE AND STABILITY

Store the test card at 4–30°C with a valid period of 24 months. Use the test card for Getein1100/Getein1160/Getein1180/Getein208 within 1 hour once the foil pouch is opened.

For test card of Getein1200/Getein1600: if the cartridge is opened, it could be stable within 24 hours once exposed to air. If the test cards can't be used up at a time, please put the

cartridge back to the foil pouch and resealed along the entire edge of zip-seal. The remaining test cards should be used up within 7 days.

PRECAUTIONS

1. For *in vitro* diagnostic use only.
2. Do not use the kit beyond the expiration date.
3. Do not use the test card if the foil pouch or the cartridge is damaged.
4. Do not open pouches or the cartridge until ready to perform the test.
5. Do not reuse the test card.
6. Do not reuse the pipet.
7. Handle all specimens as potentially infectious. Proper handling and disposal methods should be followed in accordance with local regulations.
8. Carefully read and follow user manual to ensure proper test performance.

SPECIMEN COLLECTION AND PREPARATION

1. This test can be used for **plasma and whole blood samples**. **Sodium citrate** can be used as the anticoagulant for plasma and whole blood. Samples should be free of hemolysis.
2. Suggest using plasma for better results.
3. If testing is delayed, plasma sample may be stored up to 3 days at 2–8°C or stored at -20°C for 1 month before testing (whole blood sample may be stored up to 3 days at 2–8°C).
4. Refrigerated or frozen sample should reach room temperature and be homogeneous before testing. Avoid multiple freeze-thaw cycles.
5. Do not use heat-inactivated samples.
6. **SAMPLE VOLUME (for Getein1100/Getein1160/Getein1180): 100 µL.**
(for Getein208): 60 µL.

TEST PROCEDURE

1. Collect specimens according to user manual.
 2. Test card, sample and reagent should be brought to room temperature before testing.
- For Getein1100:
3. Confirm SD card lot No. in accordance with test kit lot No.. Perform "SD card" calibration when necessary.
 4. Remove the test card from the sealed pouch immediately before use. Label the test card with patient or control identification.
 5. Put the test card on a clean table, horizontally placed.
 6. Using sample transfer pipette, deliver **100 µL** of sample into one tube of sample diluent, mix gently and thoroughly. Then drop **100 µL** of sample mixture into the sample port on the test card.
 7. **Reaction time: 10 minutes.** Insert the test card into Getein1100 and press "ENT" button or click on "Start" icon

(for Android Getein1100) after reaction time is elapsed. The result will be shown on the screen and printed automatically.

For Getein1160/Getein1180:

8. Confirm SD card lot No.in accordance with test kit lot No.. Perform "SD card" calibration when necessary.
9. Enter testing interface of Getein1160/Getein1180.
10. Remove the test card from the sealed pouch immediately before use.Label the test card with patient or control identification.
11. Put the test card on a clean table,horizontally placed.
12. Using sample transfer pipette,deliver **100 µL** of sample into one tube of sample diluent,mix gently and thoroughly.Then drop **100 µL** of sample mixture into the sample port on the test card.
13. **Reaction time: 10 minutes.** Insert the test card into Getein1160/Getein1180 immediately after sample loading. The analyzer will count down the reaction time and automatically test the card after reaction time is elapsed. The result will be shown on the screen and printed automatically.

For Getein208:

14. Long press the Power Button to start the analyzer
 15. The system will enter (Test) menu.
 16. Insert the MEMO memory chip which is with the same batch number as the test card.
 17. Select (Test) menu, press (OK) to enter [Read Calibration Card] interface.
 18. Press (OK) to automatically obtain the test item, batch number, serial number and sampling volume. Select the sample type by pressing < or > buttons.
 19. Press (OK). The screen then prompts [Insert test card] and starts counting down from 60 sec. Insert test card within the 60 sec.
- Note:** Do not move the test card after it is inserted.
20. Add sample within 120 sec when the screen prompts [Wait for sample]. Then draw **60 µL** of sample and drop it into 1000 µL of sample diluent. Then drop **60 µL** of sample mixture into the sample port on the test card.
 21. After sample adding, the system starts react-time countdown automatically.
 22. After the countdown is over, the system starts testing automatically.

Please check and record test results then.

- Note:** Test results are saved automatically in the system.
23. Long Press (OK) to return to the main interface. Take out and discard the test card.

For Getein1200/Getein1600:

24. Each cartridge for Getein1200/Getein1600 contains a specific RFID card which can calibrate automatically.
25. Place the sample diluent at the correct position in Getein1200/Getein1600.
26. Place samples in the designed area of the sample holder, insert the holder and select the right test item, Getein1200/Getein1600 will do the testing and print the result automatically.

Notes:

1. It is required to perform "SD card" calibration when using a new batch of kits for Getein1100/Getein1160/Getein1180/Getein208.
2. It is suggested to calibrate once for one batch of kits for Getein1100/Getein1160/Getein1180/Getein208.
3. Make sure the test card and the sample insertion is correct and complete.

TEST RESULTS

Getein1100/Getein1160/Getein1180/Getein208/Getein1200/Getein1600 can scan the test card automatically and display the result on the screen. For additional information, please refer to the user manual of Getein1100/Getein1160/Getein1180/Getein208/Getein1200/Getein1600.

EXPECTED VALUE

The expected normal value for D-Dimer was determined by testing samples from 500 apparently healthy individuals. The 95th percentile of the concentration for D-Dimer is 0.50 mg/L. (The probability that value of a normal person below 0.50 mg/L is 95%.)

It is recommended that each laboratory establish its own expected values for the population it serves.

PERFORMANCE CHARACTERISTICS

Measuring Range	0.10~10.00 mg/L
Lower Detection Limit	≤0.10 mg/L
Within-Run Precision	≤10%
Between-Run Precision	≤15%

LIMITATIONS

1. As with all diagnostic tests, a definitive clinical diagnosis should not be made based on the result of a single test. The test results should be interpreted considering all other test results and clinical information such as clinical signs and symptoms.
2. Samples containing interferents such as rheumatoid factor, human anti-mouse antibody and heterophile antibody may influence the results. In this case, results of this test should be used in conjunction with clinical findings and other tests. The table below listed the maximum allowance of these potential interferents.

Interferent	Hemoglobin	Triglyceride	Bilirubin
Concentration (Max)	5 g/L	25 g/L	0.1 g/L

REFERENCES






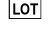







1. Sarig G, Kilil-Drori AJ, Chap-Marshak D, Brenner B, Drugan A. Activation of coagulation in amniotic fluid during normal human pregnancy. Thromb Res. 2011 Apr 18.
2. Roldán V, Marín F, Muiña B, Torregrosa JM, Hernández-Romero D, Valdés M, Vicente V, Lip GY. Plasma

von Willebrand Factor Levels Are an Independent Risk Factor for Adverse Events Including Mortality and Major Bleeding in Anticoagulated Atrial Fibrillation Patients. J Am Coll Cardiol. 2011 Apr 11.

3. Sakamoto K, Yamamoto Y, Okamoto H, Okabe M. D-dimer is helpful for differentiating acute aortic dissection and acute pulmonary embolism from acute myocardial infarction. Hellenic J Cardiol. 2011 Mar-Apr; 52(2):123-127.
4. EN ISO 18113-1:2011 In vitro diagnostic medical devices - Information supplied by the manufacturer (labelling) - Part 1: Terms, definitions and general requirements.
5. EN ISO 18113-2:2011 In vitro diagnostic medical devices - Information supplied by the manufacturer (labelling) - Part 2: In vitro diagnostic reagents for professional use.

DESCRIPTION OF SYMBOLS USED

The following graphical symbols used in or found on D-Dimer Fast Test Kit (Immunofluorescence Assay) are the most common ones appearing on medical devices and their packaging. They are explained in more details in the European Standard EN ISO 15223-1:2021.

Key to symbols used			
	Manufacturer		Use-by date
	Do not re-use		Date of manufacture
	Consult <i>instructions for use</i> or consult <i>electronic instructions for use</i>		Batch code
	Temperature limit		<i>In vitro</i> diagnostic medical device
	Contains sufficient for <n> tests		Authorized representative in the European Community/European Union
	CE mark		Do not use if package is damaged and consult <i>instructions for use</i>
	Catalogue number		

Thank you for purchasing D-Dimer Fast Test Kit (Immunofluorescence Assay). Please read this user manual carefully before operating to ensure proper use.

Version: WIF05-S-13



Getein Biotech, Inc.
 Add.: No.9 Bofu Road, Luhe District, Nanjing, 211505, China
 Tel: +86-25-68568508
 Fax: +86-25-68568500
 E-mail: tech@getein.com.cn
 overseas@getein.com.cn
 Website: www.getein.com



Lotus NL B.V.
 Add.: Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.
 E-mail: peter@lotusnl.com
 Tel: +31644168999