



Cardiac Troponin I Fast Test Kit (Immunofluorescence Assay)

IF1001 for Getein1100
IF2001 for Getein1600
IF5001 for Getein1160
IF3001 for Getein1180
IF4001 for Getein1200
IF6001 for Getein208

REF

Instructions for Use

INTENDED USE

Cardiac Troponin I Fast Test Kit (Immunofluorescence Assay) is intended for *in vitro* quantitative determination of Cardiac Troponin I (cTnI) in human serum, plasma or whole blood samples. This test is used as an aid in the diagnosis of myocardial injury such as Acute Myocardial Infarction (AMI), Unstable Angina, Acute Myocarditis and Acute Coronary Syndrome (ACS).

SUMMARY

Troponin, a molecular complex that is bound to the thin filament (actin) of striated muscle fibers, acts with intracellular calcium to control the interaction of the thin filament with the thick filament (myosin), thus regulating muscle contraction. Troponin consists of three regulatory proteins: T, which connects the troponin complex and tropomyosin (another cardiac muscle regulatory protein); I, which prevents muscle contraction in the absence of calcium; C, which binds calcium. Cardiac Troponin I (MW 22.5 kDa) and the two skeletal muscle isoforms of Troponin I have considerable amino acid sequence homology, but cTnI contains an additional N-terminal sequence and is highly specific for myocardium.

Clinical studies have demonstrated the release of cTnI into the blood stream within hours following acute myocardial infarction (AMI) or ischemic damage. Elevated levels of cTnI are detectable in blood within 4 to 6 hours after the onset of chest pain, reaching peak concentrations in approximately 8 to 28 hours, and remain elevated for 3 to 10 days following AMI. Due to the high myocardial specificity and the long duration of elevation, cTnI has become an important marker in the diagnosis and evaluation of patients suspected of having an AMI.

The current guideline of The Joint European Society of Cardiology/American College of Cardiology Committee support the use of cTnI as a preferred marker of myocardial injury. Several major studies have shown that cTnI is also a predictor of cardiac risk in patients with unstable angina. The American College of Cardiology and the American Heart Association's current guidelines recommend using troponin results when making treatment decisions regarding unstable angina and

non-ST segment elevation MI (NSTEMI).

PRINCIPLE

The test uses an anti-human cTnI monoclonal antibody conjugated with fluorescence latex and another anti-human cTnI monoclonal antibody coated on the test line. After the sample has been applied to the test strip, the fluorescence latex-labelled anti-human cTnI monoclonal antibody binds with the cTnI in sample and forms a marked antigen-antibody complex. This complex moves to the test card detection zone by capillary action. Then marked antigen-antibody complex is captured on the test line by the anti-human cTnI monoclonal antibody. The fluorescence intensity of the test line increases in proportion to the amount of cTnI in sample. Then insert test card into Getein1100/Getein1160/Getein1180 Immunofluorescence Quantitative Analyzer/Getein208 Hand-held Integrated System/automatically inserted by Getein1200/Getein1600 Immunofluorescence Quantitative Analyzer (hereinafter referred to as Getein1100, Getein1160, Getein1180, Getein208, Getein1200 and Getein1600), the concentration of cTnI in sample will be measured and displayed on the screen. The value will be stored in Getein1100/Getein1160/Getein1180/Getein208/Getein1200/Getein1600 and available for downloading. The result can be easily transmitted to the laboratory or hospital information system.

CONTENTS

1. A kit for Getein1100/Getein1160/Getein1180 contains:

- 1) Package specifications: 25 tests/kit, 10 tests/kit
- 1) cTnI test card in a sealed pouch with desiccant
- 2) Disposable pipet
- 3) User manual: 1 piece/kit
- 4) SD card: 1 piece/kit
- 5) Whole blood buffer: 1 bottle/kit

2. A kit for Getein208 contains:

- 1) Package specifications: 25 tests/kit, 10 tests/kit
- 1) cTnI test card in a sealed pouch with desiccant
- 2) Disposable pipet
- 3) Sample diluent
- 4) User manual: 1 piece/kit
- 5) SD card: 1 piece/kit

3. A kit for Getein1200/Getein1600 contains:

- 1) Package specifications: 2x24 tests/kit, 2x48 tests/kit
- 1) Sealed cartridge with 24/48 Getein cTnI test cards
- 2) User manual: 1 piece/kit

Materials required for Getein1200/Getein1600:

- 1) Sample diluent: 1 bottle/kit
- 2) Box with pipette tips: 96 tips/kit
- 3) Mixing plate: 1 piece/kit
4. Sample diluent/Whole blood buffer composition:

Phosphate buffered saline, proteins, detergent, preservative, stabilizer.

5. A test card consists of:

A plastic shell and a reagent strip which is composed of a sample pad, nitrocellulose membrane (one end of the membrane is coated with a fluorescence latex-labelled anti-human cTnI monoclonal antibody, the test line is coated with another anti-human cTnI monoclonal antibody and the control line is coated with rabbit anti-mouse IgG antibody), absorbent paper and liner.

Note: Do not mix or interchange different batches of kits.

APPLICABLE DEVICE

Getein1100 Immunofluorescence Quantitative Analyzer
Getein1180 Immunofluorescence Quantitative Analyzer
Getein1160 Immunofluorescence Quantitative Analyzer
Getein208 Hand-held Integrated System
Getein1200 Immunofluorescence Quantitative Analyzer
Getein1600 Immunofluorescence Quantitative Analyzer

STORAGE AND STABILITY

Store the test kit at 4~30°C with a valid period of 24 months. Use the test card for Getein1100/Getein1160/Getein1180/Getein208 within 1 hour once the foil pouch is opened. For test card of Getein1200/Getein1600: if the cartridge is opened, it could be stable within 24 hours once exposed to air. If the test cards can't be used up at a time, please put the cartridge back to the foil pouch and reseal along the entire edge of zip-seal. The remaining test cards should be used up within 7 days.

PRECAUTIONS

1. For *in vitro* diagnostic use only.
2. Do not use the kit beyond the expiration date.
3. Do not use the test card if the foil pouch or the cartridge is damaged.
4. Do not open pouches or the cartridge until ready to perform the test.
5. Do not reuse the test card.
6. Do not reuse the pipet.
7. Handle all specimens as potentially infectious. Proper handling and disposal methods should be followed in accordance with local regulations.
8. Carefully read and follow user manual to ensure proper test performance.

SPECIMEN COLLECTION AND PREPARATION

1. This test can be used for **serum, plasma and whole blood samples**. Heparin and EDTA should be used as the anticoagulant for plasma and whole blood. Samples should be free of hemolysis.
2. Suggest using serum or plasma for better results.
3. Serum or plasma can be used directly. For whole blood sample, one drop of whole blood buffer must be added before

testing.

4. If testing is delayed, serum and plasma samples may be stored up to 7 days at 2~8°C or stored at -20°C for 6 months before testing (whole blood sample may be stored up to 3 days at 2~8°C).
5. Refrigerated or frozen sample should reach room temperature and be homogeneous before testing. Avoid multiple freeze-thaw cycles.
6. Do not use heat-inactivated samples.
7. **SAMPLE VOLUME (for Getein1100/Getein1160/Getein1180): 100 μ L.**

(for Getein208): 70 μ L.

TEST PROCEDURE

1. Collect specimens according to user manual.
2. Test card, sample and reagent should be brought to room temperature before testing.

For Getein1100:

1. Confirm SD card lot No. in accordance with test kit lot No.. Perform "SD card" calibration when necessary.
2. Remove the test card from the sealed pouch immediately before use. Label the test card with patient or control identification.
3. Put the test card on a clean table, horizontally placed.
4. Using sample transfer pipette, deliver **100 μ L** of sample into the sample well on the test card (for whole blood sample, one drop of whole blood buffer must be added after loading **100 μ L** sample on the test card).
5. **Reaction time: 10 minutes.** Insert the test card into Getein1100 and press "ENT" button or click on "Start" icon (for Android Getein 1100) after reaction time is elapsed. The result will be shown on the screen and printed automatically.

For Getein1160/Getein1180:

1. Confirm SD card lot No. in accordance with test kit lot No.. Perform "SD card" calibration when necessary.
2. Enter testing interface of Getein1160/Getein1180.
3. Remove the test card from the sealed pouch immediately before use. Label the test card with patient or control identification.
4. Put the test card on a clean table, horizontally placed.
5. Using sample transfer pipette, deliver **100 μ L** of sample into the sample well on the test card (for whole blood sample, one drop of whole blood buffer must be added after loading **100 μ L** sample on the test card).
6. Insert the test card into Getein1160/Getein1180 immediately after sample loading. The analyzer will count down the reaction time (10 minutes) and automatically test the card after reaction time is elapsed. The result will be shown on the screen and printed automatically.

For Getein208:

1. Long press the Power Button to start the analyzer.
2. The system will enter (Test) menu.
3. Confirm SD card lot No. in accordance with test kit lot No..

- Read the relevant information in the SD card for calibration.
- Insert test card according to the analyzer prompts.
- Note:** Do not move the test card after it is inserted.
- Add sample according to the analyzer prompts. Then draw **70 µL** of sample and drop it into sample diluent. Then drop **70 µL** of sample mixture into the sample port on the test card.
 - After sample adding, the system starts react-time countdown automatically.
 - After the countdown is over, the result will be shown on the screen.

For Getein1200/Getein1600:

- Each cartridge for Getein1200/Getein1600 contains a specific RFID card which can calibrate automatically.
- Place the sample diluent at the correct position in Getein1200/Getein1600.
- Place samples in the designed area of the sample holder, insert the holder and select the right test item, Getein1200/Getein1600 will do the testing and print the result automatically.

Notes:

- It is required to perform “SD card” calibration when using a new batch of kits for Getein1100/Getein1160/Getein1180/Getein208.
- It is suggested to calibrate once for one batch of kits for Getein1100/Getein1160/Getein1180/Getein208.
- Make sure the test card and the sample insertion is correct and complete.

TEST RESULTS

Getein1100/Getein1160/Getein1180/Getein208/Getein1200/Getein1600 can scan the test card automatically and display the result on the screen. For additional information, please refer to the user manual of Getein1100/Getein1160/Getein1180/Getein208/Getein1200/Getein1600.

EXPECTED VALUE

The expected normal value for Troponin I was determined by testing samples from 500 apparently healthy individuals. The 99th percentile of the concentration for cTnl is 0.10 ng/ml. (The probability that value of a normal person below 0.10 ng/ml is 99%.)

It is recommended that each laboratory establish its own expected values for the population it serves.

PERFORMANCE CHARACTERISTICS

Measuring Range	0.10~50.00 ng/ml
Lower Detection Limit	≤ 0.10 ng/ml
Within-Run Precision	≤10%
Between-Run Precision	≤15%

LIMITATIONS

- As with all diagnostic tests, a definitive clinical diagnosis should not be made based on the result of a single test. The test results should be interpreted considering all other test results and clinical information such as clinical signs and symptoms.
- Interferents in samples may influence the results. The table below listed the maximum allowance of these potential interferents.






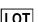

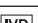





Interferent	Hemoglobin	Triglyceride	Bilirubin
Concentration (Max)	5 g/L	10 g/L	0.2 g/L

REFERENCES

- Mauro Pantaghini. Undefined International Federation of Clinical Chemistry and Laboratory Medicine (IFCC). Scientific Division Committee on Standardization of Markers of Cardiac Damage. Clin Chem Lab Med, 1998, 36:887-893.
- Antman EM, Anbe DT, Armstrong PW, et al. ACC/AHA guidelines for the management of patients with ST-elevation myocardial infarction: a report of the American College of Cardiology/American Heart Association Task Force on Practice Guidelines (Committee to Revise the 1999 Guidelines for the Manage 2004).
- EN ISO 18113-1:2011 In vitro diagnostic medical devices - Information supplied by the manufacturer (labelling) - Part 1: Terms, definitions and general requirements.
- EN ISO 18113-2:2011 In vitro diagnostic medical devices - Information supplied by the manufacturer (labelling) - Part 2: In vitro diagnostic reagents for professional use.

DESCRIPTION OF SYMBOLS USED

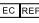
The following graphical symbols used in or found on Cardiac Troponin I Fast Test Kit (Immunofluorescence Assay) are the most common ones appearing on medical devices and their packaging. They are explained in more details in the European Standard EN ISO 15223-1:2021.

Key to symbols used			
	Manufacturer		Use-by date
	Do not re-use		Date of manufacture
	Consult <i>instructions for use</i> or consult <i>electronic instructions for use</i>		Batch code
	Temperature limit		<i>In vitro</i> diagnostic medical device
	Contains sufficient for <n> tests		Authorized representative in the European Community/European Union
	CE mark		Do not use if package is damaged and consult <i>instructions for use</i>
	Catalogue number		

Thank you for purchasing Cardiac Troponin I Fast Test Kit (Immunofluorescence Assay). Please read this user manual carefully before operating to ensure proper use.

Version: WIF02-S1-05

 Getein Biotech, Inc.
 Add.: No.9 Bofu Road, Luhe District, Nanjing, 211505, China
 Tel: +86-25-68568508
 Fax: +86-25-68568500
 E-mail: tech@getein.com.cn
 overseas@getein.com.cn
 Website: www.getein.com

 CMC Medical Devices & Drugs S.L.
 Add.: C/ Horacio Lengo N° 18, CP 29006, Málaga, Spain
 Tel: +34951214054